

Open size : 364mm x 495mm
Closed Size : 85mm (w) x 85mm (l) x 320mm (h)
Material : ArtCard 300 gsm with Matt Lamination



PENILE PUMP
INSTRUCTIONS

How To Use

- 1. Ensure all parts are securely connected.
- 2. Apply a small amount of lubricant to the silicone/TPE ring.
- 3. Gently insert your penis, making sure the seal is airtight.
- 4. Press the “+” button to turn on the pump.
 - Each press increases the pressure range.
 - The pump automatically cycles pressure in a “pulse” pattern for safety

Pressure Range	Value (mbar)
Level 1	–50 to –125
Level 2	–125 to –200
Level 3	–200 to –275
Level 4	–275 to –350

- 5. To release pressure, press the “valve” button.

Before Use

- 1. Wash the pump with warm water and antibacterial soap to disinfect.
- 2. Use a **water-based lubricant** (e.g., KY Jelly).
- 3. Do **not** use massage oils, hand creams, detergents, or chemicals as lubricants.

After Use

- Clean immediately with warm water and antibacterial soap.
- Air dry naturally, avoid direct sunlight.
- Store in a cool, dry, dust-free place.
- Keep out of reach of children (small parts may pose a choking hazard).

For your best result, please use this pump with a sex lubricant. Always select the appropriate pressure to avoid hurting your genital.



USB
Rechargeable



Multi-Speed
Control



Clinical-Grade
Technology



Pro-Level
Intensity




Silent Step



Wet-Dry Ready



EMPOWER
YOUR
PERFORMANCE

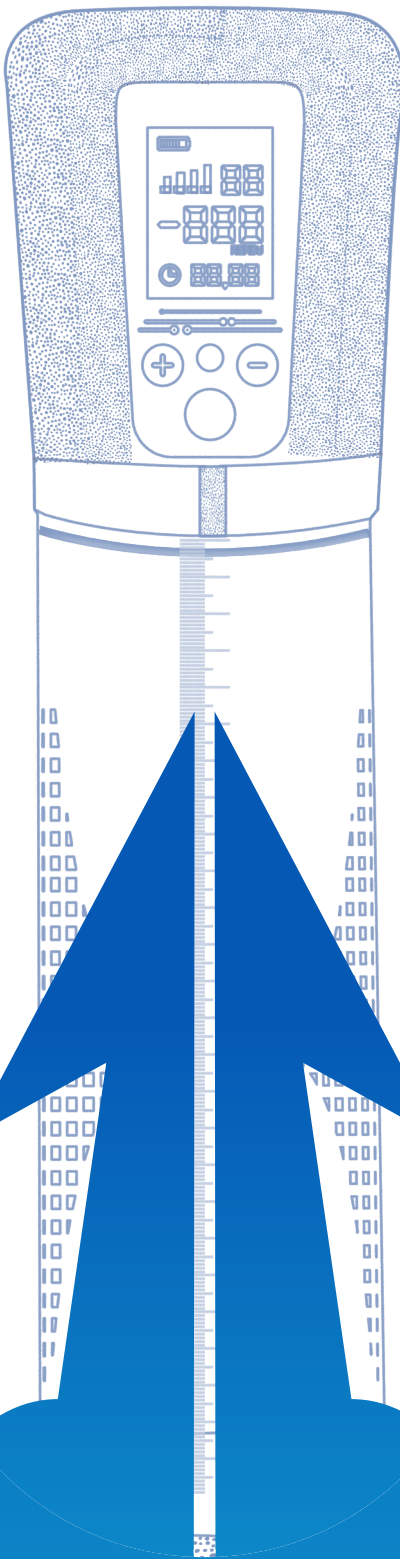

MENSS redefines men’s

wellness with precision

engineering, medical-grade

safety, and discreet design

Engineered for performance.

Designed for confidence.